

Integra® Miltex®

Sterilization Containers

Limit uncertainty by protecting
your instrument investment



INTEGRA®
LIMIT UNCERTAINTY

Summary of Integra Sterilization Container Product Features

Gaskets within retention plates and container perimeter ensure a secure seal of the container, lid and filter.

Snap on filter plate with a simple locking mechanism for quick, secure and easy change of disposable filters.

Round disposable filters with centering hole to assist in proper alignment.

Lid and Bottom manufactured from anodized aluminum alloy. Easy to maintain, corrosion resistant, lightweight yet durable, and conducts heat for material drying period.

4 ID Tag Windows (2 front & 2 back)

ID tags available in 6 colors and can be engraved with up to 18 characters per line.


Low profile latch easily and securely closes lid over container bottom.

Tamper Evident Lock – identifies sterility breach.

Fold-down handles extend and stop at 90° for a secure grasp to aid lifting and avoid pinching.

Corner design allows for same size container stacking in storage.

2 Indicator Card Windows (1 front & 1 back). Indicator card has visible bar for clear indication that steam or gas sterilization has been achieved.

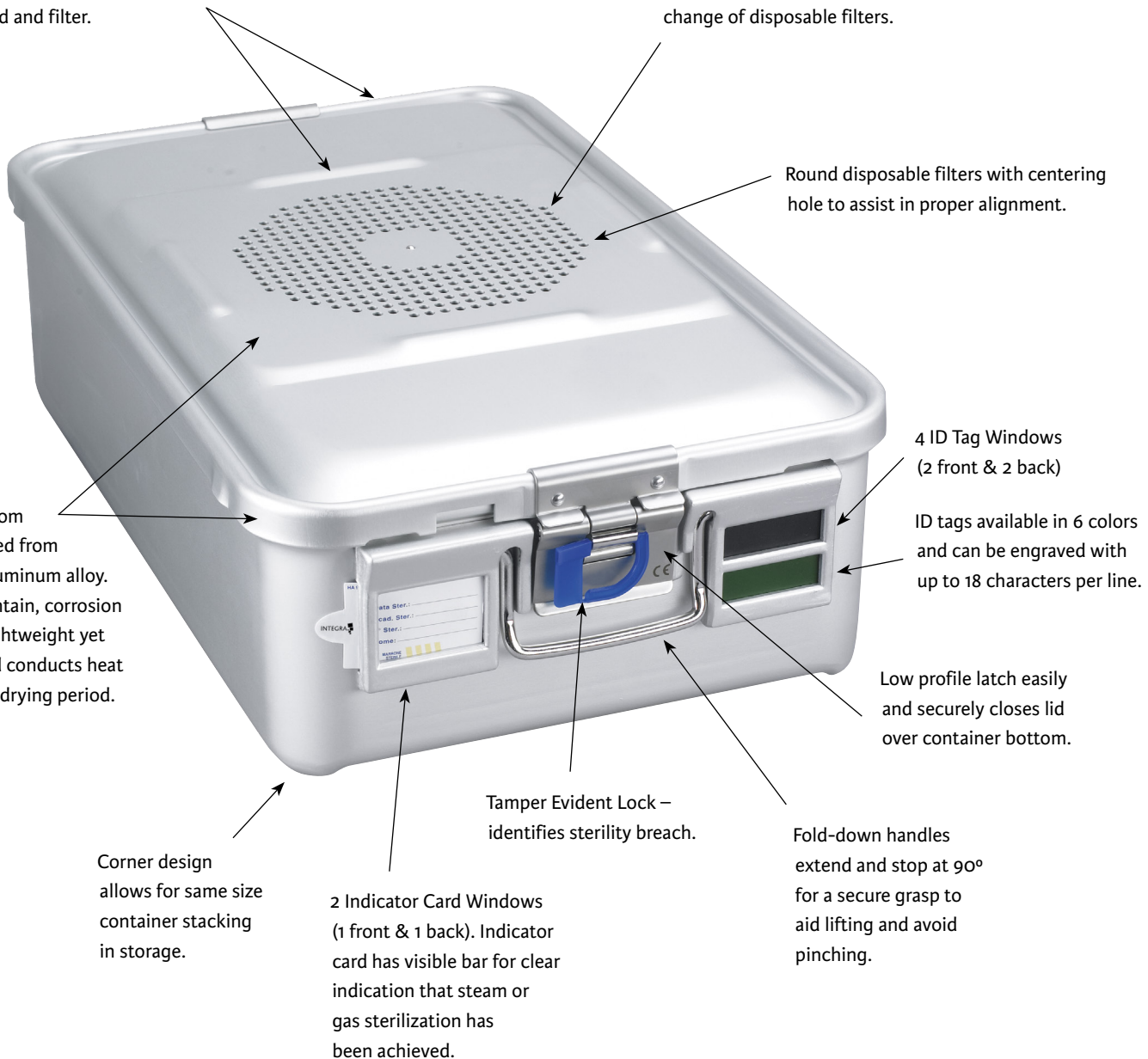

*Baskets are constructed with solid wire mesh to allow for maximum ultrasonic cleaning; handles can be positioned to assist with a secure grasp for removal under sterile conditions.

One-Half Size (11¼" x 11") Containers

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom

Height	Solid Bottom Full Container	Solid Bottom Only	Perforated Bottom Full Container	Perforated Bottom Only
4"	730400	730400B	730410	730410B
5½"	730401	730401B	730411	730411B
6"	730402	730402B	730412	730412B
8"	730403	730403B	730413	730413B
10"	730404	730404B	730414	730414B

Note: Complete container assembly only available in silver.
See below for other color options with lids and bottoms sold separately.



One-Half Size Colored Lids (bottoms sold separately)

Includes Filter Retention Plate and Latches

Color	Lid Only
Silver	3610000
Blue	3610001
Red	3610002
Gold	3610003
Green	3610004
Black	3610005



One-Half Size Wire or Perforated Baskets

Height	Wire Basket	Perforated Basket
2"	740400	750400
3"	740401	750401
4"	740402	750402
Lid		750420



Wire Basket



Perforated Basket

Three-Quarter Size (18¼" x 11") Containers

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom

Height	Solid Bottom Full Container	Solid Bottom Only	Perforated Bottom Full Container	Perforated Bottom Only
4"	730500	730500B	730510	730510B
5½"	730501	730501B	730511	730511B
6"	730502	730502B	730512	730512B



Three-Quarter Size Colored Lids (bottoms sold separately)

Includes Filter Retention Plate and Latches

Color	Lid Only
Silver	3630000
Blue	3630001
Red	3630002
Gold	3630003
Green	3630004
Black	3630005




Three-Quarter Size Wire or Perforated Baskets

Height	Wire Basket	Perforated Basket
2"	740500	750500
3"	740501	750501
4"	740502	750502
Lid		750520



Wire Basket

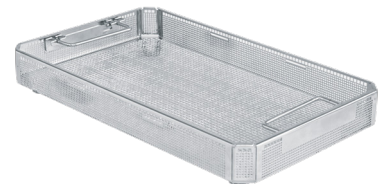

Perforated Basket

Full-Size (23" x 11") Containers

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom

Height	Solid Bottom Full Container	Solid Bottom Only	Perforated Bottom Full Container	Perforated Bottom Only
4"	730600	730600B	730610	730610B
5 1/3"	730601	730601B	730611	730611B
6"	730602	730602B	730612	730612B
8"	730603	730603B	730613	730613B
10"	730604	730604B	730614	730614B

Note: Complete container assembly only available in silver.
See below for other color options with lids and bottoms sold separately.



Full-Size Colored Lids (bottoms sold separately)

Includes Filter Retention Plate and Latches

Color	Lid Only
Silver	3650000
Blue	3650001
Red	3650002
Gold	3650003
Green	3650004
Black	3650005




Full-Size Wire or Perforated Baskets

Height	Wire Basket	Perforated Basket
2"	740600	750600
3"	740601	750601
4"	740602	750602
Lid		750620



Wire Basket

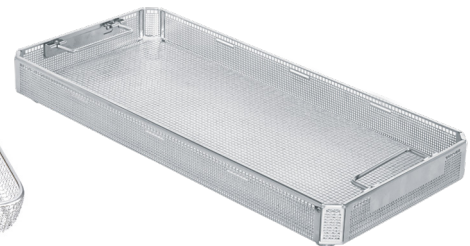

Perforated Basket

Wide-Body (23¹/₃" x 15³/₄") Containers

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom

Height	Solid Bottom Full Container	Solid Bottom Only	Perforated Bottom Full Container	Perforated Bottom Only
120mm	730701	730701B	730711	730711B
180mm	730702	730702B	730712	730712B



Wide-Body Lid (bottoms sold separately)


Includes Filter Retention Plate and Latches

Color	Lid Only
Silver	730220



Wide-Body Wire Baskets

Height	Wire Basket
2"	740700
4"	740701

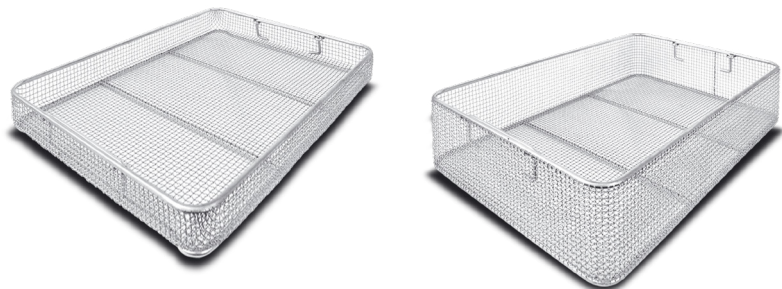

Identification Tags

Each tag allows for 2 lines of custom text, 18 characters per line.

Catalog Number	Color
775100	Silver
775101	Blue
775102	Red
775103	Gold
775104	Green
775105	Black




Snap-on Basket Tags (with key ring)

Each small tag allows for 3 lines of custom text, 32 characters per line.

Each large tag allows for 5 lines of custom text, 32 characters per line.

Color	Small	Size	Large	Size
Silver	775120	2 3/8" x 1/2"	775130	2 3/8" x 1 1/4"
Gold	775121	2 3/8" x 1/2"	775131	2 3/8" x 1 1/4"
Red	775122	2 3/8" x 1/2"	775132	2 3/8" x 1 1/4"
Blue	775123	2 3/8" x 1/2"	775133	2 3/8" x 1 1/4"
Green	775124	2 3/8" x 1/2"	775134	2 3/8" x 1 1/4"
Black	775125	2 3/8" x 1/2"	775135	2 3/8" x 1 1/4"

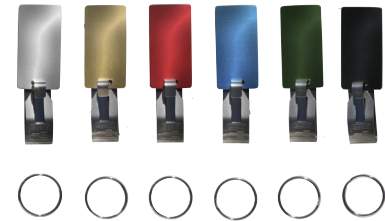


Tamper Evident Locks and Indicator Cards

Catalog Number	Description
775400	Tamper Evident Locks, single use, blue, box of 1000
775300	Indicator Cards, box of 250




Tamper Evident Lock

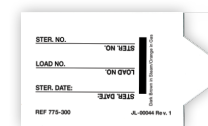


Indicator Card

Retention Plate and Filters

Catalog Number	Description
775212	Retention Plate
775200	Disposable paper filters, single use, 7½", round, box of 1000, GRAVITY
775250	Disposable paper filters, single use, 7½", round, box of 1000, PRE-VACUUM

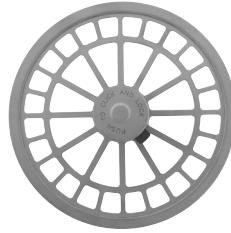

Retention Plate

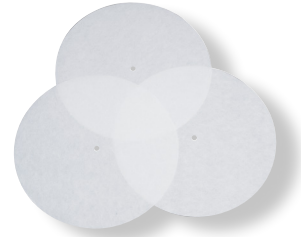

Disposable Paper Filters

Silicone Holders

Catalog Number	Description	Size
775630	Silicone Holder #1	8mm x 26mm x 40mm
775631	Silicone Holder #2	8mm x 26mm x 40mm
775632	Silicone Holder #3	8mm x 26mm x 40mm
775633	Silicone Holder #4	10mm
775634	Silicone Holder #5	5mm
775650	Silicone Holder Clamps, 4 slot, plastic	



#1



#2



#3




#4



#5



Silicone Mats

Catalog Number	Description
776400	OneHalf Size Silicone Mat
776500	ThreeQuarter Size Silicone Mat
776600	FullSize Silicone Mat

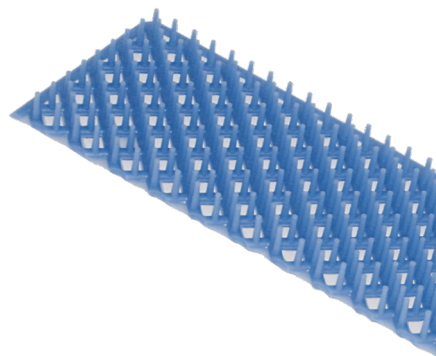

Partition Sheets and Holding Pins

Catalog Number	Description	Size
775610	Partition Sheet	2" x 7/8" (50mm x 20mm)
775611	Partition Sheet	5" x 7/8" (130mm x 20mm)
775612	Partition Sheet	9 7/8" x 7/8" (225mm x 20mm)
775613	Partition Sheet	18 1/8" x 7/8" (460mm x 20mm)
775600	Holding Pin for Wire Baskets	
775601	Holding Pin for Partition Sheets	

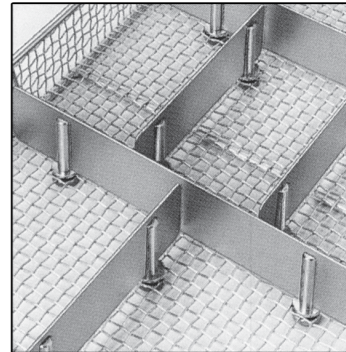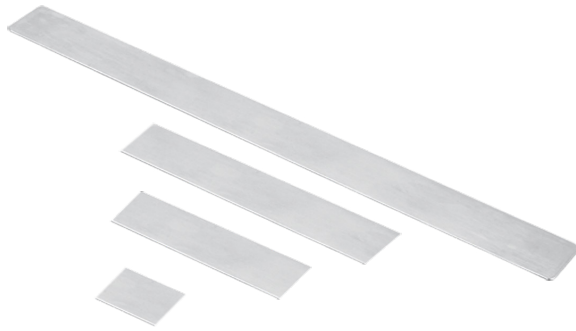

Holding Pin for Wire Baskets



Holding Pin for Partition Sheets

Holding Clamps

Catalog Number	Size Range
775620	1/8" to 1/2" (7mm-12mm)
775621	1/2" to 5/8" (12mm-16mm)
775622	5/8" to 1" (16mm-26mm)
775623	1" to 1 1/2" (26mm-36mm)
775624	1 1/2" to 1 3/4" (36mm-45mm)



Tolerances in product specifications, due to the hand craftsmanship of the instruments, are subject to minor variations in dimensions and finishes.

For more information or to place an order, please contact:

Integra ▪ 589 Davies Drive, York, PA 17402
877-468-5572 USA ▪ 914-789-7094 outside USA
609-750-4259 fax ▪ 914-592-8056 outside USA fax
integralife.com

